

KaVo QUATTROcare®

KaVo QUATTROcare®. Perfect care.



KaVo. Dental Excellence.

Easy to use.

Given the daily tasks of my practice, there is little time left to care for my high quality KaVo instruments. There is not much time left for the office team either. It is therefore, all the more important to rely on a well-conceived care system that anyone can operate. Insert the instruments. Close the front door. Start. Finished. Four at once. Just like that. Excellent.

Just a press of the button. Fully automated. The KaVo QUATTROcare® offers maximum speed and user-friendliness. The KaVo QUATTROcare® is so easy to operate. Simultaneously insert four instruments on the four couplings, close the door and start the care program. A complete care cycle lasts only 60 seconds. The KaVo QUATTROcare® is immediately ready to use.



- 1 *Insert the instruments.*
To insert the instruments, a MULTIflex coupling is available as well as three INTRAmatic couplings.
- 2 *Close the door.*
And you're ready to go.
- 3 *Press the start button*
Fully automated cleaning and care of up to four instruments in just 60 seconds.
- 4 *Finished.*

Optimum care. Maximum efficiency.

The difference is in the details. KaVo QUATTROcare® cleans even the smallest niche. For safe removal of protein (RKI recommendation in Germany - prevention of infections in dentistry), prior cleaning with KaVo CLEANspray and DRYspray is required. The speed of the KaVo QUATTROcare® lessens the workload in my practice and gives me the assurance of perfect care. The service life of my instruments is extended.

KaVo QUATTROcare® cleans by foaming the cleaning components. In addition, the automated dosing system ensures optimum lubrication of the turbines, handpieces and contra-angle handpieces.

The KaVo QUATTROcare® was developed based on the most current quality and safety standards. And it is universal. That means, you can clean your high-quality KaVo instruments using the KaVo QUATTROcare® as well as instruments from other manufacturers.

A long service life is a question of cleaning and care and not just quality. With the KaVo QUATTROcare®, KaVo has set another international benchmark for efficiency in instrument care. The innovative KaVo QUATTROcare® system that has been proven for years is incomparably thorough, economical and easy to use. Perfect for all my instruments that are used daily in my practice.

"Dental Excellence" led you to select KaVo turbines, handpieces and contra-angle handpieces. The KaVo QUATTROcare® is the logical extension of this ideology. To your great benefit. The internationally established KaVo care system is optimised for the daily requirements

of your practice, and joins the latest care technology with highly modern design and maximum efficiency. KaVo QUATTROcare® maintains the unique quality of the instruments and extends their life regulations.



- 1 *The KaVo QUATTROcare® effect.*
The forceful expansion of the QUATTROcare® cleaning agent inside the instruments dissolves and removes residue.
- 2 *Efficiency in the dental practice.*
KaVo SUPERtorque LUX 650 B, cleaned with the KaVo QUATTROcare®. Seven years in daily use without a breakdown or exchanging the rotor. Always optimally clean and maintained, always ready to use.



It's in your hand.

Every KaVo instrument is top quality. And it will stay that way.

To keep my KaVo instruments from losing their precision and value, I care for them with the KaVo QUATTROcare®. At the same time, I create the ideal preconditions for additional hygiene. You can only save costs and achieve perfect results with perfectly cleaned and hygienically prepared instruments.


The KaVo QUATTROcare® comes in two designs. As a standard unit, or with the innovative purging function. Compressed air additionally blows through the spray channels to remove any remaining oil residue and thoroughly clean the KaVo instruments.

Moreover, the device is available in two colours: mint-green or anthracite.

Advantages at a glance:

- 1 minute care cycle for 4 instruments
- Program-controlled automatic care cycle
- Perfect interior cleaning from both the aerosol foaming/ expansion effect and of course the aerosol's propellant pressure
- Optimum lubrication with automated dosing
- Easy operation with the push of a button prevents errors
- The spray channels are blown clean with the purging function*
- Elimination of excess oil usage, prevents costly wastage and unnecessary mess
- Basic setup includes 4 couplings (1 x MULTIflex/3 x INTRAmatic)
- Universally applicable (Adapters for all conventional makes)
- Environmentally-friendly KaVo QUATTROcare® plus spray
- Modern and ergonomic design
- Available in mint-green or anthracite

*only available with the KaVo QUATTROcare® 2104 A

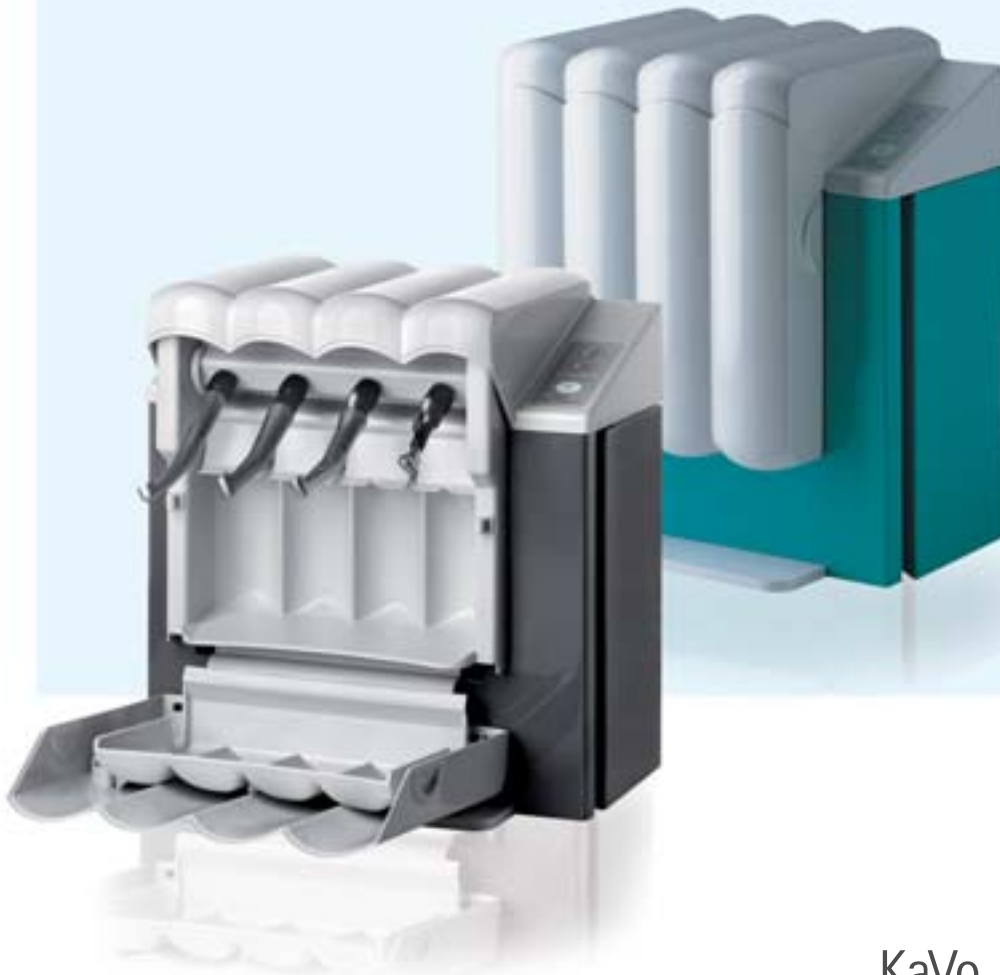

*KaVo QUATTROcare® 2104.
Device in the standard design without
compressed air connection*

*KaVo QUATTROcare® 2104 A.
Device with compressed air connection to
blow clean the smallest spray channel*



KaVo. Dental Excellence.