

KaVo. Dental Excellence.



GENTLEsilence – low noise with high torque.



- > Comfortable for both dentist and patient thanks to very low noise level during treatment (57 dB(A))
- > Best visual access thanks to very small head dimensions (H = 13.1 mm; \emptyset = 12.5 mm)
- > More space thanks to the patented KaVo 100° / 19° head / knee angle combination
- > Save working time by highly efficient preparations with high performance (8000 B: 19 W; 8000 BS: 17 W; 8000 BN: 18 W)
- > Improved hygiene: the integrated anti-reflux system prevents ingress of contaminated water spray

8000 B Mat.-Nr. 1.000.0800 8000 C Mat.-Nr. 1.000.2800 8000 BS Mat.-Nr. 1.000.5600 8000 BN Mat.-Nr. 1.000.5700



Quality, Longevity, Relia & Value-for-Money

GENTLEsilence 8000 B/C: Double protection for a long service life



KaVo *Master ^{Series}*

A true master knows no compromise.

- Highest treatment quality
- Maximum treatment comfort
- Highest patient satisfaction



KaVo **Expert** ^{Series}

Real experts know their strengths.

Best price-endurance ratio

- Proven configuration solutions
- Highest efficiency




KaVo. Dental Excellence.

KaVo GENTLEsilence 8000 B/C

Endurance is all in the head.

KaVo GENTLEsilence 8000 B/C now with steel head casing and anti-reflux system.


KaVo **Master** Series

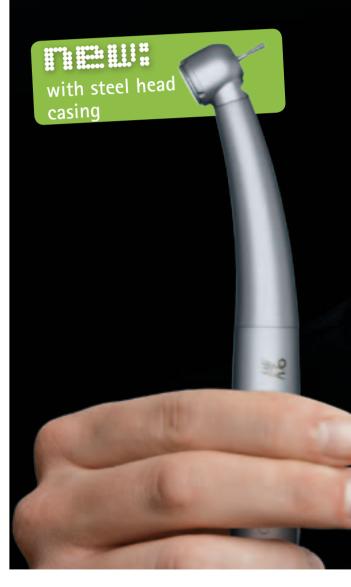

Even more durable thanks to the steel head casing.

- > Head casing made from high-quality steel
- > Impervious to corrosion
- > Robust against mechanical damage

GENTLEsilence 8000 with steel head casing – extremely durable by comparison.

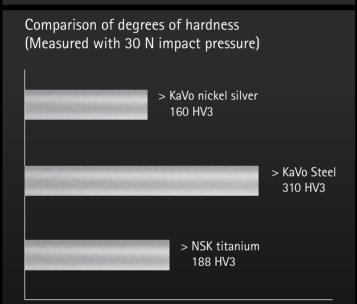

Degree of hardness



Maximum durability due to double protection.


1. protection: steel head casing

> GENTLEsilence 8000 with steel head casing – for steely long-term protection.



2. protection: anti-reflux system

 > GENTLEsilence with anti-reflux system – for hygienic long-term protection.

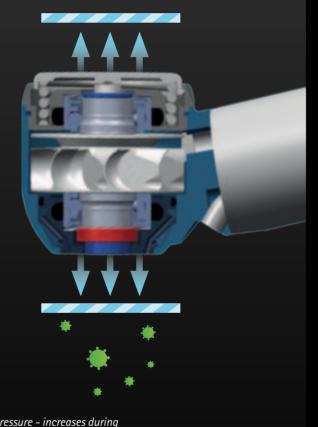


Overpressure - increases during slowing down phase of the cartridge

Even more durable thanks to the

anti-reflux system.

- > Optimum conservation of value as there is little ingress of contamination into the interior of the handpiece
- > Best specification for effective reprocessing and improved hygiene

GENTLEsilence 8000 B – minimal aspiration properties.

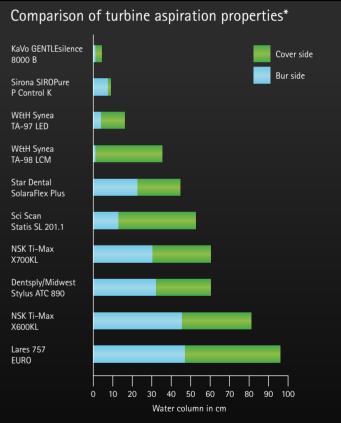

* Sealed measuring under vacuum during the slowing down phase (measured in cm of water)



KaVo. Dental Excellence.