


KaVo GENTLEmini LUX 5000 B

Compact and Powerful –
The KaVo GENTLEmini is the
ideal addition to your
instrument set.

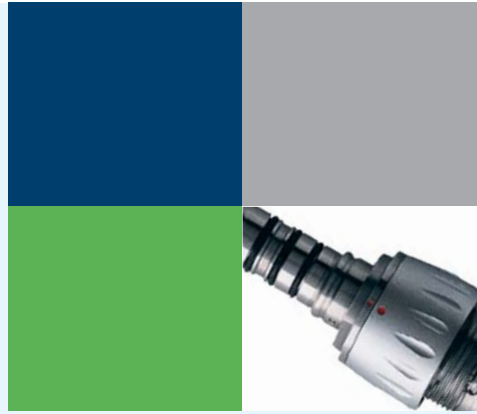

KaVo. Dental Excellence.

The smaller head height allows for superior range of visibility.

The new KaVo Turbine with miniature head is designed specifically to allow improved access especially when treating children and seniors

The GENTLEmini is ideal when preparation conditions are restricted principally in the molar region or when treating patients with smaller mouths. The miniature head offers optimised visibility and hence the right conditions for unhindered work even when space is restricted. The GENTLEmini LUX turbine with its smaller head dimensions can substantially improve treatment results. The small head leaves more room to work and optimum access and is the perfect addition to your existing instrument set.



GENTLEmini LUX 5000 B.


*with 16 mm miniature drill bit.




GENTLEsilence LUX 8000 B.

Compact and powerful – Main features and advantages for you:

- Head height, head diameter and patented 100° head angle – Everything is optimized for superior visibility and unhindered work, even when access is restricted.
- 17 W output offers a highly effective, powerful performance without compromise.
- Maximum security that you can rely on – the KaVo chucking system with a strong 29 N grip.
- As strong and gentle as the GENTLEsilence:
Quiet, vibration-free operation at only 59dB (A)) for pleasant and relaxed work.
- Quality that pays – the KaVo original ceramic ball bearings ensure a long service life due to proven technology.

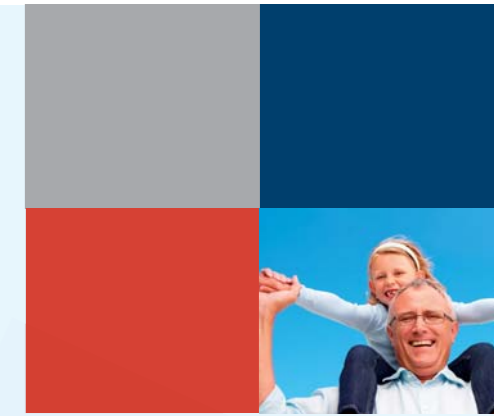

Just one click: The proven KaVo MULTIflex coupling is the basis for the KaVo GENTLEmini LUX. The KaVo GENTLEmini LUX can be used on any treatment unit.

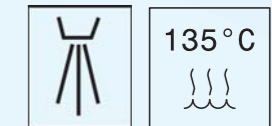

The GENTLEmini LUX can be washed and disinfected in a thermodisinfector and sterilized up to 135°C in an autoclave, and thereby satisfies the latest hygiene requirements.

The three-nozzle spray and exchangeable micro-filter ensure optimum and reliable cooling.

The reduced head height of the GENTLEmini LUX can be optimally used with 16 mm miniature drill bits. 19 mm FG drill bits can also be used without a problem.

The precise illumination of the treatment area with 25,000 lux ensures optimum visibility.

The Plasmatec coating used by KaVo provides an excellent grip and outstanding hygiene.

Automatic drive pressure regulation ensures continuous, ideal operating conditions and a long service life.



KaVo GENTLEmini LUX 5000 B
Order number: 1.002.3300



Lasting value from optimum care quality.

KaVo handpieces have a longstanding reputation for long service life.

KaVo instruments reflect superior quality. And it will stay that way. To ensure that the KaVo instruments do not lose any of their attractiveness and value, we offer a perfect care program tailored to KaVo products.

Superior quality is ensured, and the traditionally long service life is maintained.

From KaVo QUATTROcare to KaVo SPRAYrotor to KaVo Spray.

For high-quality KaVo instruments, KaVo also offers high-quality care products.

These are optimally tailored to the required care of the individual instruments, and the STERlclave is an ideal supplement for maintaining all of the current hygiene requirements.

KaVo SPRAYrotor 3 N

Cleans and maintains your valuable instrument set in a highly-effective, economical and environmentally friendly manner with automatic dosing.



KaVo Spray

For manual care and cleaning of all KaVo handpieces and contra-angle handpieces, heads and turbines.

STERlclave

With a cycle duration of 15 - 36 minutes (controlled programs), the KaVo STERlclave is significantly faster than conventional sterilisers.



KaVo QUATTROcare 2104 A

Lasting value par excellence. The innovative care system automatically cleans and maintains up to four instruments simultaneously.



KaVo. Dental Excellence.