

KaVo MULTIflex® System

An extremely versatile connection.

KaVo MULTIflex® is not only a coupling but also the key to the entire system. Because the KaVo MULTIflex® coupling is the connection to an extensive range of performance perfect KaVo instruments for all indications.

KaVo air motor 181 M and hence all KaVo shanks, heads and straight handpieces.

KaVo INTRAfLEX® LUX motor instruments

KaVo SONICflex® airscaler, the oscillating instrument for a diverse and incredibly far reaching range of clinical applications.

All KaVo turbines with MULTIflex® connection.

More than you think can be mounted on the MULTIflex® coupling:

KaVo PROPHYflex® 3 for professional tooth cleaning.

KaVo CORONAfLEX® for easy removal of dental prostheses.

KaVo RONDOflex®, the superior alternative to conventional air abrasion instruments.

KaVo QUATTROcare® KaVo SPRAYrotor® KAVOspray®

Systematic care.

Care can be a complicated matter. However, regular care is a decisive factor in the life of any dental instrument – irrespective of the quality features of the instrument in terms of material and manufacture.

KaVo provides you with suitable care options, which are effective and easy to use. They are optimally tailored to meet the precise needs of KaVo product care.


KaVo SPRAYrotor 3N Intensive cleaning with Powered rotation into the smallest handpiece recess.

KaVo QUATTROcare 2104/A Innovative care system, automatically cleans and cares for four of your valuable instruments at the same time. Fast as a flash.

KAVOspray 500 ml of efficiency. Manual care for all contra-angle shanks, handpieces, heads and turbines from KaVo.

KaVo POWERtorque turbines 646 B/C.

Just a click that's the trick.

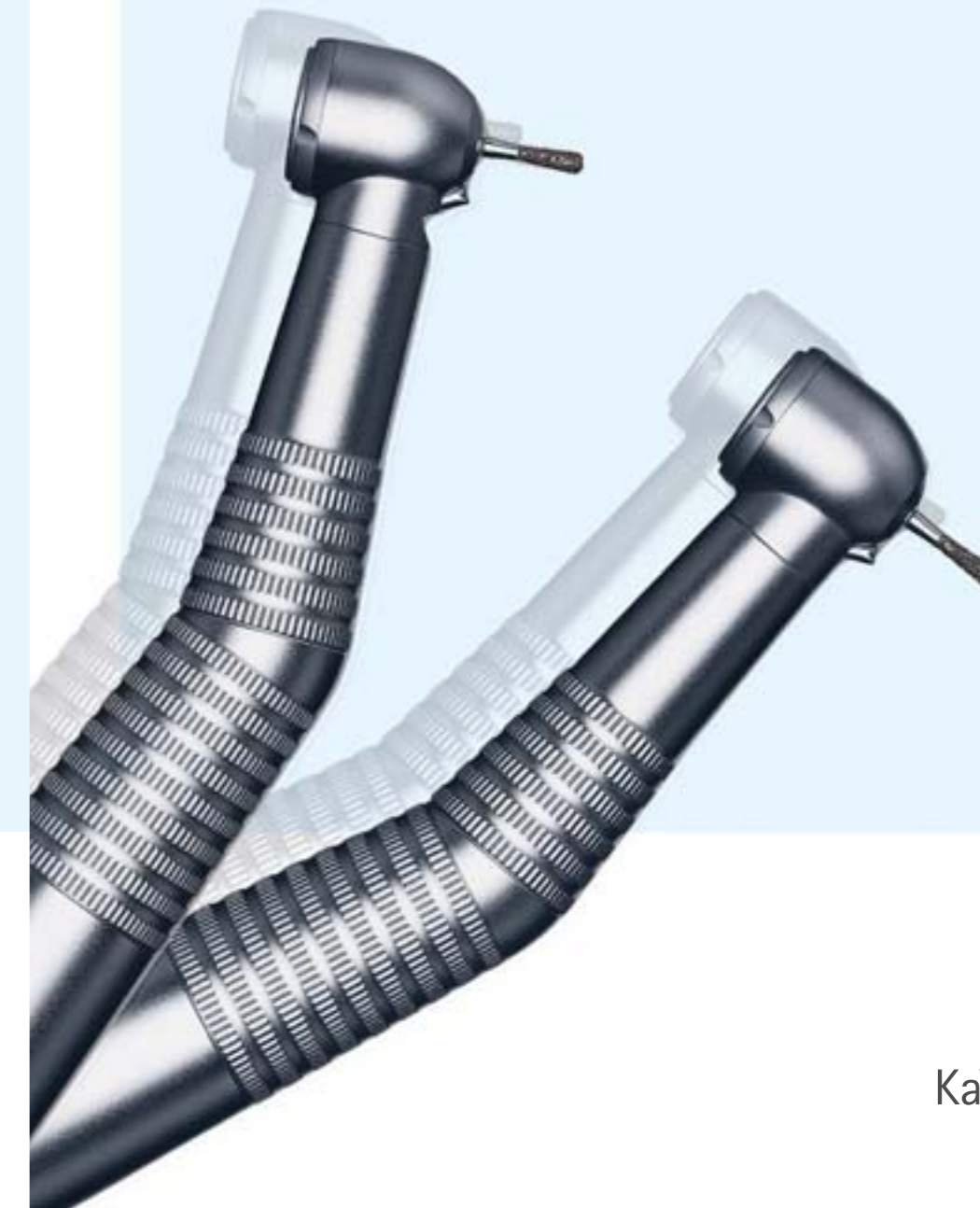

Mat.-Nr. 1.004.12.0 11/03-11/10 We reserve the right to make technical modifications. Slight colour differences are due to the printing process. © Copyright KaVo Dental GmbH.



KaVo. Dental Excellence.

KaVo Dental GmbH · D 88400 Biberach/Riß · Telefon +49 7351 56-0 · Fax +49 7351 56-1488 · www.kavo.com



KaVo. Dental Excellence.


MULTiflex® LUX coupling 460 LE. MULTiflex® coupling 460 E.

The extra-fast connection.

The KaVo MULTiflex® coupling saves you time-consuming interruptions for instrument change during your work. Because one action - click! - is sufficient to fix the new KaVo instrument on the MULTiflex® coupling. Without a great deal of fiddling. Without loss of time.

Once fitted the handpiece is rotatable through 360 degrees. There is no troublesome tension in the wrist and no twisting of the tubing. The built-in KaVo suck-back stop prevents ingress of contaminated spray water. A good link and an extra-fast connection for your successful practice!

KaVo MULTiflex® LUX 460 LE
REF: 1.001.7599

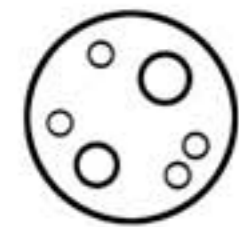

4-hole standard connection with light

KaVo MULTiflex® 460 E
REF: 1.001.7600



4-hole standard connection without light

KaVo POWERtorque LUX 646 B. KaVo POWERtorque 646 C.

The new standard for your practice.

The KaVo POWERtorque turbine will set new standards in your practice. Because POWERtorque combines various advantages that optimally support your work.

The turbine is powerful and at the same time quiet. It removes a great deal of material but still runs completely smoothly. It has a small head for compactness in the oral cavity.

It has proven precision ball bearings and can be perfectly disinfected and sterilized. It feels great externally and rests securely in the hand.

You can decide between POWERtorque and the POWERtorque LUX, which has the 3rd generation fibre optic eye.



KaVo original precision ball bearings.
The proven technology for safe work and a long cartridge life. The vibration-free running permits precise work with little pressure.




Always cool.
The single port fan-shaped spray ensures constant optimal cooling of the handpiece and of the entire preparation area. The spray wets down to the bur tip and is particularly finely distributed.




POWERtorque 646 C without light.
REF: 1.001.7523


POWERtorque LUX 646 B with light.
REF: 1.001.7522

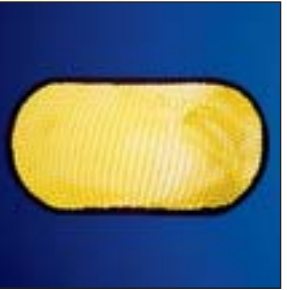

Truly illuminating.
The illuminance of 25,000 lux of the POWERtorque LUX leaves nothing in the dark. The high-precision strong glass rod ensures excellent light transmission and a very bright illumination of the preparation area even after many thousands of sterilisation cycles and disinfection applications.

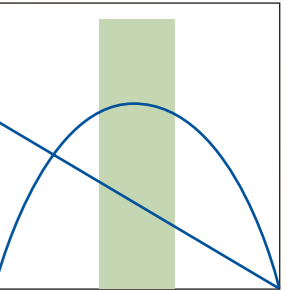

Safe preparation.
High power of 15 watt - that's power with plenty in reserve even for difficult preparations with removal of copious quantities of tooth tissues and other material. Power is decisive.

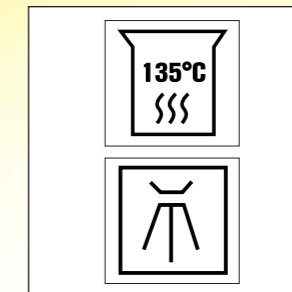

Easy to care for.
Sterilization and a long life are not mutually exclusive. All KaVo turbines can be sterilized at 135°C in an autoclave and washed externally in a Miele thermal disinfectant.



Completely quiet.
Precision requires a steady hand. And vibration-free running. At every speed range.



Very safe.
The proven KaVo pushbutton chuck system ensures secure tension and easy changing of the bur.